

A UNIQUE DEVICE
FOR HIGH QUALITY
MEASUREMENT,
IN A
TOUCH



M'EYE FIT TOUCH™

DIGITAL MEASURING DEVICE

M'eye Fit™
TOUCH

A UNIQUE TOUCH IN MEASUREMENT



PROFESSIONAL

It brings valuable benefits straight to the sales table, including quality measurement and electronic content. M'eye Fit Touch is a genuine professional device designed for demanding opticians, equipped with a robust hybrid-PC architecture and a high-resolution camera.

SIMPLE

No need to raise the camera at eye level, or take photos from multiple angles. With the customer seated at your desk, you achieve precise measurement in two shots. It's simple, fast, and compatible with all types of frames.



PRECISE

Measurement quality is comparable to all premium devices on the market today. M'eye Fit Touch uses a patented 3D measurement protocol which takes into account the exact angle of the camera and compensates for parallax. It covers PD, half PDs, fitting heights, pantoscopic angle, wrap, vertex and reading distances.

RELIABLE

On wearers, all tests show exact same value on M'eye Fit Touch and pupillometers: reproducibility is less than 0.5mm in more than 95% of the cases.* A flash secures corneal reflect.

*Result based on 600 measures on real wearers.

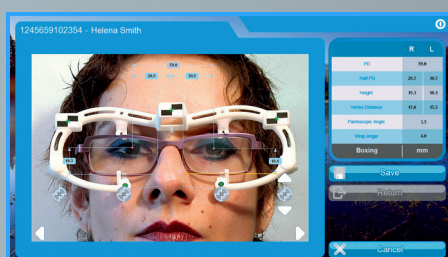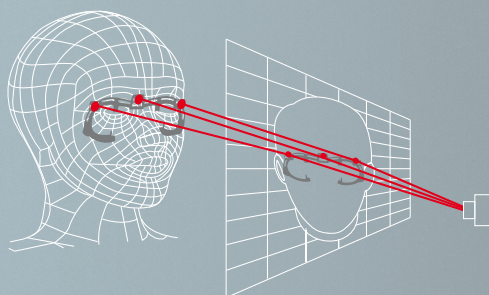

	R	L
PD	68.4	68.4
half PD	34.2	34.2
Height	69.2	69.2
Vertex Distance	42.8	42.8
Pantoscopic Angle	13.1	13.1
Wrap Angle	14.1	14.1
Bowing	mm	mm

A PERFECT FIT TO YOUR NEEDS



ENHANCED DIALOG WITH YOUR CUSTOMERS

- M'eye Fit Touch fits perfectly at your sales table like no other device.
- Seated in front of your customers, you can share measurement results easily thanks to the convenient swivel stand.
- When it's time to discuss lenses, use the educational module to show simulations and display the benefits more easily.



AN ULTRA-FLEXIBLE SALES TOOL

- M'eye Fit Touch is portable to go along with the customer in the shop.
- The integrated camera makes it easier to choose the right frames from the shelves.
- It integrates several applications.



EASY INTEGRATION

- M'eye Fit Touch is equipped with its own internal database.
- It is ready to share data with your ordering and Practice Management Systems.

ENJOY A UNIQUE NEW APPROACH TO EYE CARE!

With M'eye Fit Touch, you're fully equipped to significantly enhance the overall customer experience!

IN BRIEF

- High-quality patented measurement protocol
- Simple and user-friendly
- Easy swivelling
- Ready for data exchange

MEASUREMENTS

- PD, half PDs
- Fitting heights
- Pantoscopic angle
- Vertex distance
- Wrap*
- Reading distance

*Measurement by wrap angle template.

ACCESSORIES

- Clip
- Anti-slip stand
- Digitizer pen
- Micro SD Card + adaptator
- Power supply
- Cleaning cloth

SPECIFICATIONS

- 300(L) x 200(D) x 230(H) mm
- Weight: 1010g
- Power supply:
110/240v - 50/60Hz
- Screen:
11.6-inch touch screen
- CE conformity marking

M'eye Fit Touch is a trademark of Essilor International.

To allow ongoing improvements, these specifications are non-contractual and are subject to modification without notice.



ESSILOR INTERNATIONAL Instruments Department

64 bis, avenue Aubert – 94306 Vincennes Cedex France – Tel.: 33 (0)1 72 94 71 00 – Fax: 33 (0)1 72 94 70 27